



# **Industries and application areas**

VERDERAIR HC-PURE hygienic double diaphragm pumps are designed for hygienic applications in various industries. They can handle a wide range of low and high viscous fluids which makes them ideal for numerous applications.



## Beverage

Fruit juice, flavours, colouring.



Ingredients, alcohol, chemicals.



## Cosmetics

Shampoo, cremes, lotions, soap, toothpaste, perfume.



## Dairy industry

Milk, cream, yoghurt, cheese, butter, ice cream.



## Food

Soups, starches, sauces, pastes, jams, pie filling, dough, oil, bakery ingredients.



## Brewery/winery

Beer, wine, wort, yeast, flavouring.

## **HC-PURE** range overview


The Verderair HC-PURE hygienic double diaphragm pump is available in two designs: the FD series for hygienic applications at the beginning of the production process and the SB series for hygienic applications in the final stage of the production process.



- ✓ Superior finish:
- Avoids bacterial trapsIncreases cleanability



- ✓ Rotatable stand:
- Easy drainableQuick cleaning

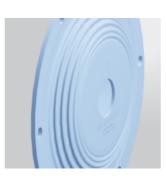

- ✓ Overmolded diaphragms:
- Easy servicingLonger lifetime
- Maintenance free air valve:
   100% self-lubricating
  - No contamination





	Material center section	FD	SB
Concept		Hygienic series	High hygienic series
Max. flow [I/min]		600	650
Max. discharge pressure [bar]		7	7
Max. operating temperature (°C)	PE PA	70 90	70 90
Max. solids particles [mm]		12	12
Internal surface roughness [µm]		< 3.2 μm	< 0.8 μm
Certifications		FDA EX	FDA EX

HYGIENIC DOUBLE DIAPHRAGM PUMPS

# VERDERAIR HC-PURE

Hygienic durability due to ultimate cleanability



- ✓ Best cleanability
- ✓ Easy to maintain
- ✓ Highest hygienic design principles





Utrechtseweg 4A
3451 GG Vleuten
The Netherlands

MAIL info@verderliquids.com
WEB www.verderliquids.com

AUSTRIA / BELGIUM / BULGARIA / CHINA / CZECH REPUBLIC / FRANCE GERMANY / HUNGARY / INDIA / ITALY / THE NETHERLANDS / POLAND REPUBLIC OF KOREA / ROMANIA / SLOVAKIA / SOUTH AFRICA SWITZERLAND / UNITED KINGDOM / UNITED STATES OF AMERICA

# Verderair HC-PURE: Hygienic durability due to ultimate cleanability

The Verderair HC-PURE series of hygienic double diaphragm pumps is designed for various hygienic applications. Their cleanability capabilities reduces the risk of bacteria traps and contamination. The pumps are manufactured according to the highest hygienic design standards for hygienic applications. Easy maintenance and long lifetime results in less downtime for your production process.

## **Best cleanability**

Verderair HC-PURE pumps are electropolished. This results in a lower average roughness which makes products more difficult to adhere and avoids bacterial traps and increases the cleanability. In addition, the surface is given greater resistance to corrosion.

The pumps can be cleaned via strip cleaning or CIP/SIP cleaning. You can easily drain the pump by turning it in the rotatable stand.



## Easy to maintain

The housing screws and manifold clamps are easily removed for complete disassembly and cleaning. The PTFE/EPDM overmolded diaphragms have a long lifetime together with the maintenance free air valve which is 100% self-lubricated results in less maintenance and downtime.

The pumps can be equipped with



## **Highest hygienic** design principles

Design and production is according to the highest hygienic design principles. HC-PURE pumps are compliant with

The single piece manifold and smooth corners and bends results in no crevices for bacteria.

hygienic DIN 11851 or Tri-clamp



EC01935/2004 and FDA.



- Leak dection
- Stroke sensor

**Options** 

→ Barrier system

## **Key features**

- 1 Electropolished SS 316L
- Solid machined center block PE/Polyamide
- 3 Rotatable stand for quick draining & service
- 4 Single piece manifolds wich smooth bends and corners results in no crevices for bacteria
- 5 PTFE/EPDM overmolded diaphragms
- 6 Manifold clamps for easy disassambly
- Maintenance free air valve

















